

ChIP Validated H3K14ac (Clone RM130) Antibody with Positive Primer Set



Catalogue no: 900022

Chromatrap®'s ChIP Validated H3K14ac Antibody and Positive Primer Set provides a complete set of tools to assist with a successful ChIP assay. Including: H3K14ac antibody, control rabbit IgG, and a positive primer set. The ChIP Validated H3K14ac Antibody and Positive Primer Set is not suitable for use with non-human species.

Background:

Histone 3 (H3) is one of the core histone proteins, comprising the protein component of chromatin. H3 is ubiquitous within chromosomes and can be found bound to most gene sequences throughout the genome. Acetylation of histone H3 occurs at several different lysine positions on the histone tail and is catalysed by a group of enzymes called histone acetyltransferases (HATs). The histone acetylation mark H3K14ac can indicate enhancer sites; helping to identify active enhancers from poised ones. H3K14ac is associated with active transcription and can be found at the start site of actively transcribed genes.

A rabbit IgG is included in this Antibody Primer Set as a negative control for the ChIP experiment.

The H3K14ac positive primer set recognises the promoter of a gene associated with active transcription and is a suitable target for this antibody.


Suggested Usage:

Component	Suggested Dilution	Figure
H3K14ac	2:1 (antibody: chromatin)	1
Rabbit IgG	2:1 (antibody: chromatin)	1
Positive Primer Set	Dilute from 4µM (provided) to 1µM working concentration	

Please note: Optimal dilutions should be determined by the user. These volumes are stated as guidelines only.

Fig 1. H3K14ac ChIP qPCR

Chromatin immunoprecipitation (ChIP) assays were performed using the Chromatrap® standard ChIP spin column sonication kit for qPCR (Cat no. 500071) with 1µg of chromatin from Hec50 cells and 2µg of Anti-H3K14ac antibody. qPCR was used to analyse the enrichment of H3K14ac onto the positive gene target locus compared with the negative gene target locus.

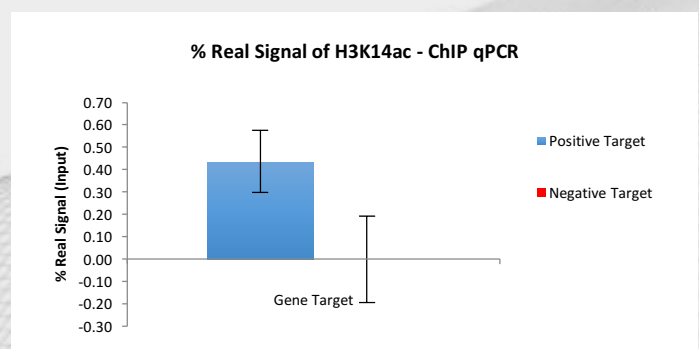

*This product is for research use only. There is a possibility that results may vary between antibody lots.



www.chromatrap.com

Clywedog Rd South
Wrexham Industrial Estate
Wrexham LL13 9XS, United Kingdom

Tel: +44 (0) 1978 666239/40
Email: support@chromatrap.com

Applications: ChIP

Concentration: 1mg/ml

Size: 50 μ l

Specificity: Human

Source: Rabbit

Type: Monoclonal

Purification: Protein A (affinity purified)

Storage Conditions: The H3K14ac antibody should be stored at +4°C. Rabbit IgG and the positive primer set should be stored at -20°C (*Avoid multiple freeze/thaw cycles as this may denature the antibody and degrade the primer sets*)

Advancements in Epigenetics

*This product is for research use only. There is a possibility that results may vary between antibody lots.

Company registered in England no 2047818 Registered office: 7 Regis Place, Bergen Way, Kings Lynn, PE30 2JN United Kingdom